



GE HealthCare

T2100-ST Series treadmills

Advanced performance and reliability
in exercise testing



T2100-ST Series treadmills

The T2100-ST Series¹ treadmills feature an extra-long, and extra-wide walking surface, as well as a zero-start feature that gradually increases speed for patient comfort and smooth operation. Ensure device and procedural accuracy—and clinical confidence.

Optimized design

- With fewer moving parts than other treadmills, the T2100-ST Series treadmills deliver smooth, quiet operation that facilitates accurate stress testing acquisition—even at high workloads
- Designed low-to-the floor, the T2100-ST Series offers a cushioned deck for patient comfort and safety
- The speed range* of 0.1 to 15.0 mph (0.2 to 24.0 km/h) self-calibrating and adjustable in 0.1 mph (0.1 km/h) increments helps to maintain protocol accuracy and repeatability even when testing heavy patients

Stopping methods

The treadmills can be stopped using several different methods to achieve the desired result:

- Decelerate slowly to a gradual stop using the emergency stop button or stop quickly using the device keypad or the tether worn on the patient's wrist
- All stop modes are available at all times giving the technician more control
- The optional adjustable height front handrail provides easy access to various size patients as needed

Easy interface

The T2100-ST Series treadmills were designed as part of the GE HealthCare stress exercise solution, and it communicates bi-directionally with GE HealthCare Stress ECG systems. This communication means both devices are constantly exchanging operating commands and data. It can easily interface with:

- CASE™ Cardiac Assessment System for Exercise
- CardioSoft™ Diagnostic System multifunction PC-based diagnostic software

Safety features

There is a prominent emergency stop button and a patient safety tether can be installed on either the left or right treadmill handle. Other safety features include:

- Dual comparative speed sensors
- Auto runaway shutdown
- Auto communication loss shutdown
- Manual twist lock Emergency Stop button
- Manual Stop Tether
- Braking system for safe patient off-loading
- Fire-rated motor pan hood enclosure

T2100-ST2/T2100-ST1 treadmills specifications

Weight	Patient weight capacity: 500 lb (227 kg) Treadmill net weight: 425 lb (193 kg)
Drive system	Heavy duty 4 hp (6 hp peak) hi-torque drive, brushless, DC servo motor
Power requirements	T2100-ST1: 100 to 130 VAC, 1-phase, 60 Hz, 20-amp power supply Must have dedicated 20 A service T2100-ST2: 200 to 264 VAC, 1-phase, 50 to 60 Hz, 15-amp power supply Must have dedicated 15 A service
Speed range*	0.1 to 15.0 mph (0.2 to 24.0 km/h) Self-calibrating and adjustable in 0.1 mph (0.1 km/h) increments
Incline range	0 to 25%, 0.5% Incremental movements, self-calibrating



Running surface	22 × 63 in (56 × 160 cm) MasterTrack® running belt tracking system Cushioned running deck absorbs shock of foot falls Self-lubricated and reversible running deck Step-up, low profile deck height: 7 in (18 cm) from floor
Chassis material	All steel with powder coat finish
Communication ports	RS232 Female Serial port USB 1.0 “B” port
Floor surface footprint	33.0 × 78.5 in (84 × 200 cm) Level surface

Operating and storage condition recommendations

Operating temperature range	+40° to +85°F (+4.5° to +38°C)
Storage temperature range	-40° to +158°F (-40° to +70°C)
Operating and storage relative humidity range	10% to 90%, Non-condensing
Altitude	-50 to 5,280 ft (-15.24 to 1,609.34 m)

Optional features:

Adjustable height handrail	31.2 to 37.5 in (79 to 92.5 cm) As measured from the top of the running surface
----------------------------	---

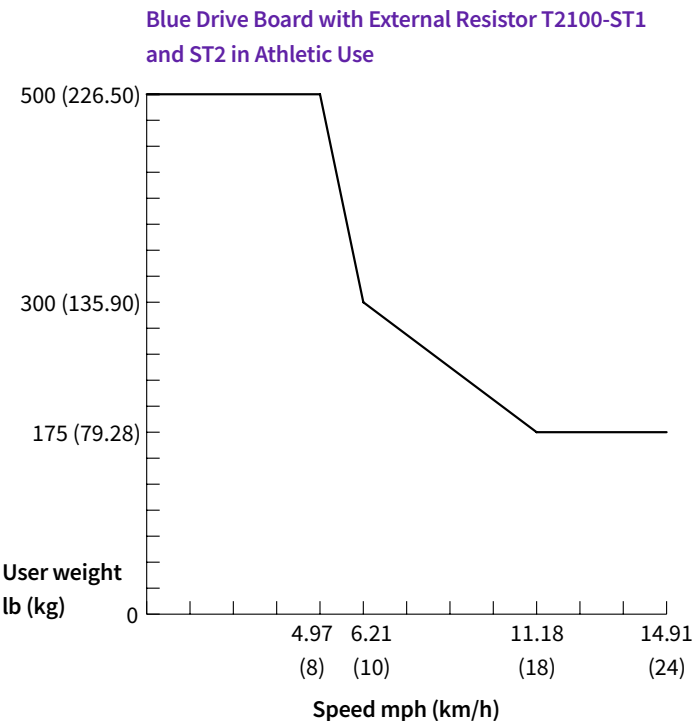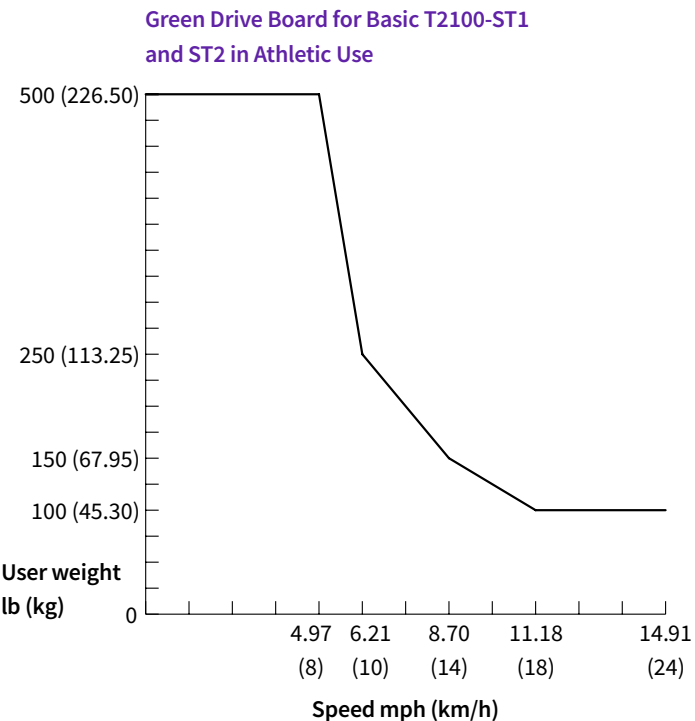

* Refer to the performance graphs.

1. The T-2100-ST is available in 2 versions, T2100-ST1, 100 to 130 VAC and T2100-ST2, 200 to 264 VAC—see the Power Requirements for full information.

© 2026 GE HealthCare

Not all products or features are available in all markets. Contact a GE Healthcare Representative for more information. Please visit www.gehealthcare.com. Data subject to change.

GE HealthCare reserves the right to make changes in specifications and features shown herein, or discontinue the product described at any time without notice or obligation. Contact your GE HealthCare Representative for the most current information. Trademarks belong to their respective owners. CASE and CardioSoft are trademarks of GE HealthCare. GE is a trademark of General Electric Company used under trademark license. MasterTrack is a trademark of Full Vision Inc.